


Updating SureTrend 4

Introduction

This article will provide step-by-step instructions on how to back up your database and install the current version of SureTrend 4.


Solution

To verify your version of SureTrend 4, check the **About** section under the **Admin** tab. The current version of SureTrend 4 is 4.0.3.4 – Build 11. If you do not have the current version, proceed with the following steps. If you have the current version, your syncing problems are related to another issue. Contact Hygiena Technical Support using by clicking [here](#).

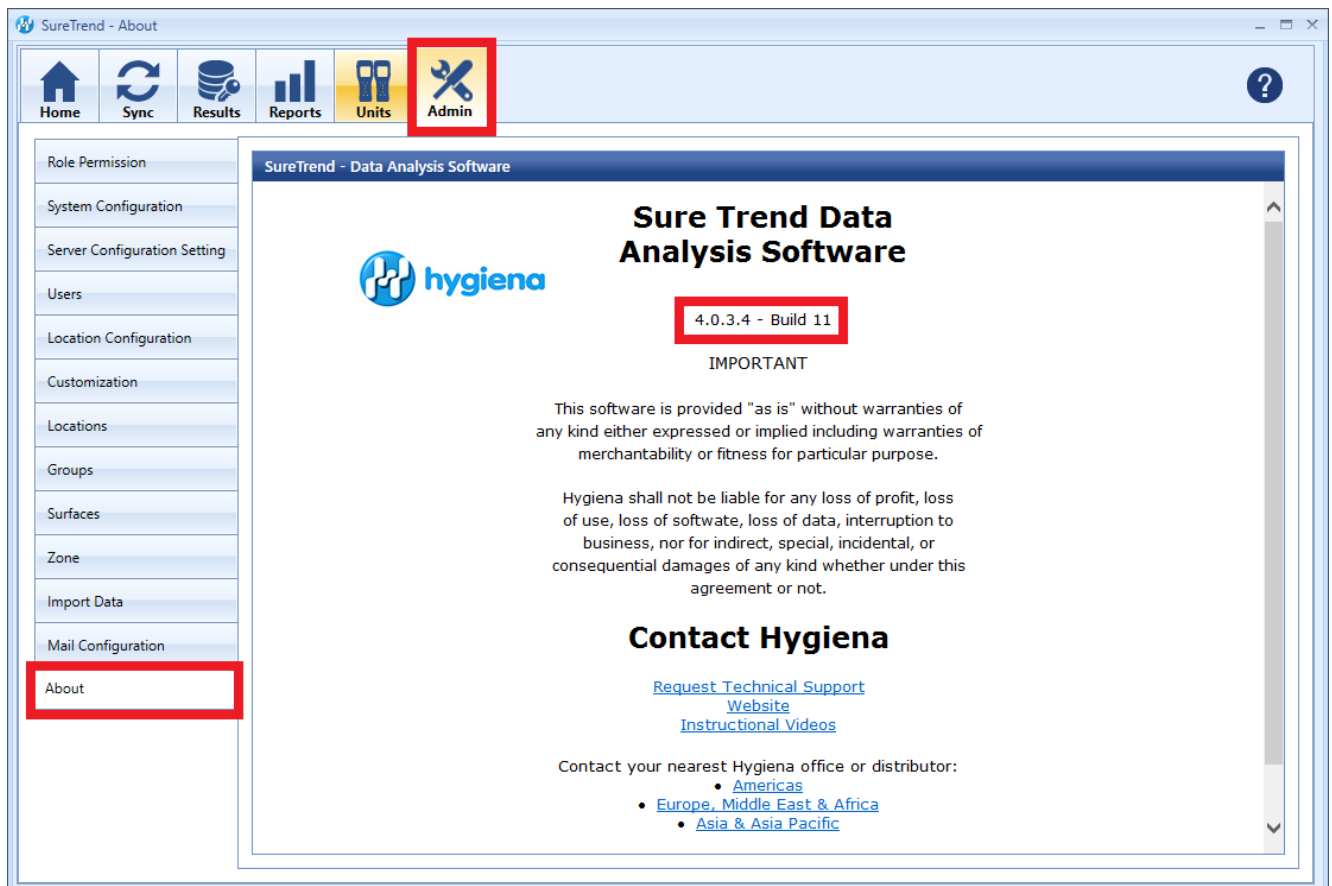


Create a Copy of Your SureTrend 4 Database

1. Navigate to your SureTrend 4 database. By default, it is located at C:\ProgramData\Hygiena\SureTrend

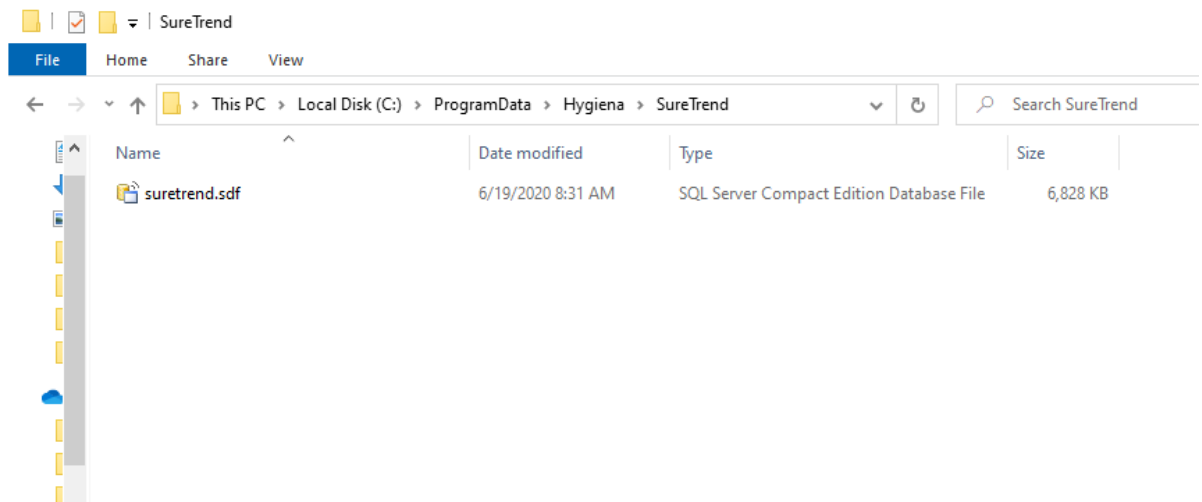

2. If you do not see your ProgramData folder, open File Explorer. Click **View** and check **Hidden Items** in the **Show/Hide** tab.

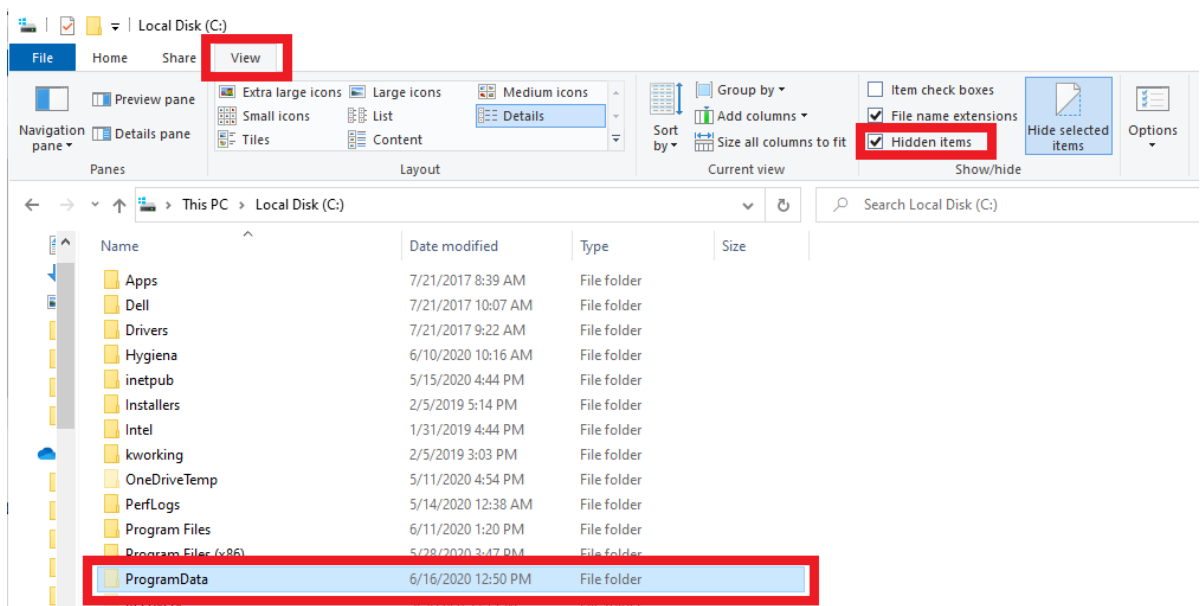

3. Copy the *.sdf database file and paste it in another location outside of this folder.

Install the Current Version of SureTrend 4

1. Download the most recent version of SureTrend 4 [here](#).
2. Run the installer, select your preferred language, and click **Next**.
3. For a standard installation leave the Multiuser Distributed Network checkbox unchecked.
4. For a Multiuser Distributed Network using internet information services

check the checkbox. (Rare)

5. Click **Finish** and launch SureTrend 4.
6. Click **Yes** when prompted to update the database.
7. You should be able to sync any new instruments to SureTrend 4 now.

Contact Hygiena Technical Support for further assistance.

- Phone: 1-888-HYGIENA (1-888-494-4362, option 2)
- Email: techsupport@hygiena.com
- [Submit a Support Ticket](#)
- [Schedule a Microsoft Teams meeting with support](#)